


E A Q R S

FOREVER YOUNG
FREE OF GLASSES IN MOST OF DAILY ACTIVITIES

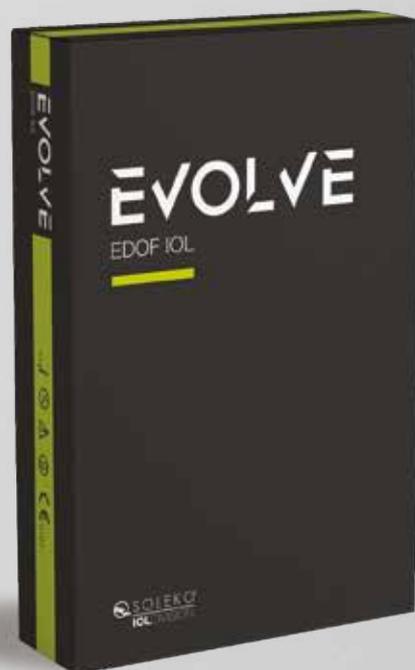

 **SOLEKO™**
IOLDMISION

EVOLVE

FREE OF GLASSES IN MOST OF DAILY ACTIVITIES

DESIGNED TO DELIVER A COMFORTABLE VISION TO EVERY PATIENT, ALSO WITH LIMITED PUPIL APERTURE, THE EVOLVE OFFERS AN INNOVATIVE EXTENDED DEPTH OF FOCUS (EDOF) FOR FAR AND INTERMEDIATE VISION.

EVOLVE EDOF

OPTIC DIAMETER

6.0 mm

TOTAL DIAMETER

11.2 mm to 11.8 mm
(depends on selected diopter)

HAPTIC ANGULATION

5°

MATERIAL

foldable acrylic with 25% H₂O and UV filter

REFRACTIVE INDEX

1,461 (546 nm , 20°C in water)

DIOPTR RANGE

from -5.00 to +30.00 (step 0.5D)

RECOMMENDED INJECTOR SYSTEM

Medicel 2.1 or equivalent

RECOMMENDED CONSTANTS

A constant (ultrasound biometry) 118.7

A constant (immersion biometry) 119.0

Optical biometry:


SRK/T IOL constant 119.1


SRKII IOL constant 119.2

Holladay I SF: 1.90

Haigis a0: 0.051, a1: 0.140, a2: 0.197

Hoffer Q/Holladay II pACD: 5.68

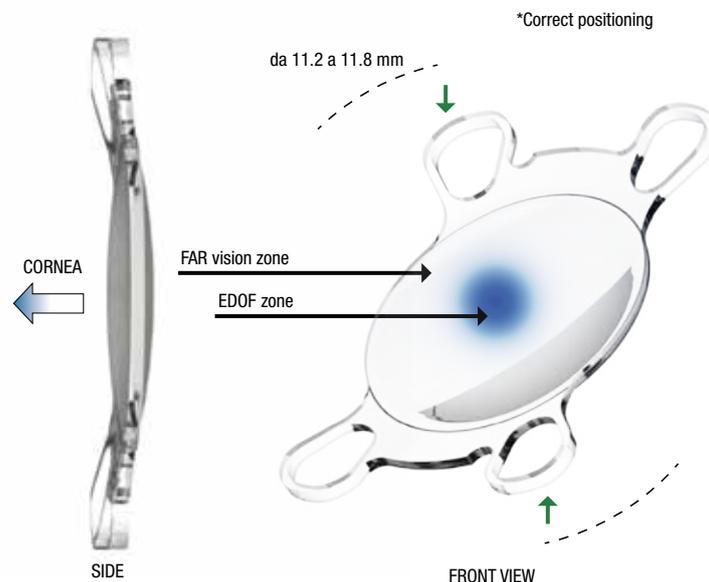

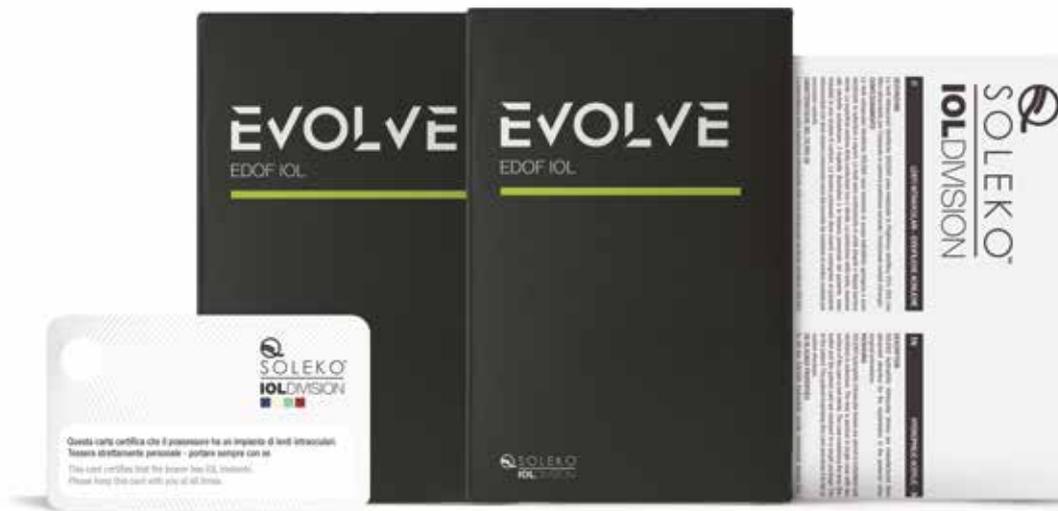

EVOLVE TORIC CUSTOMIZED

AVAILABLE DIOPTERS* (SPHERE)

STEP 0.5 D


from -5.00 to +9.50 D
from +10.00 to +19.50 D
from +20.00 to +25.00 D
from +25.50 to +35.00 D

AVAILABLE CYLINDERS* (CYL)

STEP 0.25 D AXIS 0-180, STEP 1°

from 1.00 to 10.00 D
from 1.00 to 15.00 D
from 1.00 to 10.00 D
from 1.00 to 6.00 D

*Correct positioning

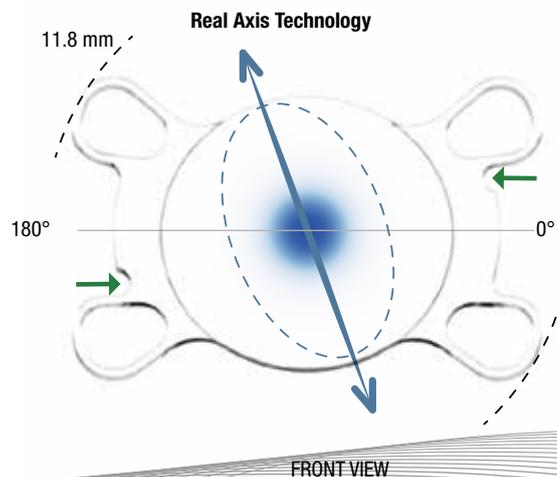

*On the lens label it is reported the spherical equivalent ($SE = sph + 1/2 cyl$). For out of range spheres or cylinders please contact ioldivision@soleko.it

Real Axis Technology
The cylinder axis is set during the construction. The lens must always be positioned at 0°-180°.

Allora cominciai la fabbricazione di alcuni oggetti necessari.
For the next ten or twelve days we sailed close to the shore.

Sommario di giornale | 20 pt.

La sera, nelle valli al centro dell'isola, ho scoperto che vi pascolano tante capre.
They signed to me to come ashore, showing they would bring us some food.

Libro di narrativa | 15 pt.

Questo pensiero produsse un particolare effetto sul mio cuore e riempi i miei occhi di lacrime.
I also found a bottle of rum which I took with me. There was no other boat on the ship

Libro tascabile | 12 pt.

Il fenomeno mi ha lasciato così sgomento, che sono rimasto in silenzio, immobile, sbalordito, come morto.
Then I convinced myself that maybe I had seen the mark of my own foot and I had been a fool to get scared by it.

Articolo quotidiano | 10 pt.



Soleko S.p.A.
03037, Via Ravano snc
Pontecorvo (FR) - ITALY
info@soleko.it
MADE IN ITALY

www.soleko-iol.it

